

Wheelchair platform scale KERN MWB

Compact wheelchair platform scale with two integrated access ramps for easy access



NEW



Your KERN specialist dealer



Features

- Especially suitable for weighing patients in wheelchairs, because the low-profile platform can be approached from either side
- The large steel weighing plate also makes it ideal for weighing obese patients
- Secure and non-slip positioning with height-adjustable rubber feet
- Level indicator to level the balance precisely
- Hold function:
When patients do not stand or sit completely still, a stable weight is calculated using an average weight and this is then "frozen". This means that you have sufficient time to attend to the patient first and then take the weight reading in peace
- BMI function to determine underweight/normal weight/surplus weight
- Wall mount for the display device
- The scale can be easily transported using the two handles and two rollers and does not require much storage space


Technical data

- Large backlit LCD display, digit height 25 mm
- Dimensions of display device
W×D×H 200×130×60 mm
- Overall dimensions,
W×D×H 1050×1050×100 mm
- Batteries included, 6×1,5V AA, operating time up to 40 h
- Cable length of display device approx. 300 cm
- Net weight approx. 55 kg

Accessories

- Protective working cover over the display device, standard, can be re-ordered, scope of delivery: 5 items, KERN MBC-A06S05
- Cleaning cloths, alcohol-free cloths for wipe disinfectant. Fulfill the legal requirements for occupational safety in accordance with TRGS 525/540. Packaging contents 80 pcs., size 20×22 cm per cloth, KERN MYC-01
- For more details see the internet

STANDARD

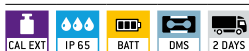

Model	Weighing range	Readout	Weighing surface	
	[Max] kg	[d] kg	mm	
KERN MWB 300K-1	300	0,1	800×800	

Pictograms



Data interface RS-232:

To connect the balance to a printer, PC or network.



Protection against dust and water splashes IPxx:

The type of protection is shown in the pictogram



Battery operation:

Ready for battery operation. The battery type is specified for each device.



Weighing principle: Strain gauges

Electrical resistor on an elastic deforming body.



Pallet shipment:

The time required for internal shipping preparations is shown in days in the pictogram.